

NEWSLETTER

Service

Issued by : PMS-DMC Hamburg
Publication No. : 4512 980 59881

Reference No. : XRD Hamburg 074
Date : May 2005 - reissue

Information for :  FSE  Managers  Specialists

Earth Leakage current circuit breaker trip off caused by a defective Varistor (SIOV)

High-speed / dual-speed rotor control for OPIMUS RAD / RF / C


Trip-off of earth leakage current breaker caused by the varistor of the rectifier filter PCB

The varistor (SIOV EY200 R13) between N and PE has a protective function.
The varistor is intended to take up any surge which might cause the destruction of the rotor control unit.

When a certain threshold value is exceeded the earth-leakage current breaker (FI switch) is released; both components (varistor and FI) are part of the safety concept.

The specification of a varistor refers to:

- Energy storage 33Joule for 2ms
- Continuous load 1W
- Surges may occur during switching operation or lightning. A surge is defined as follows:

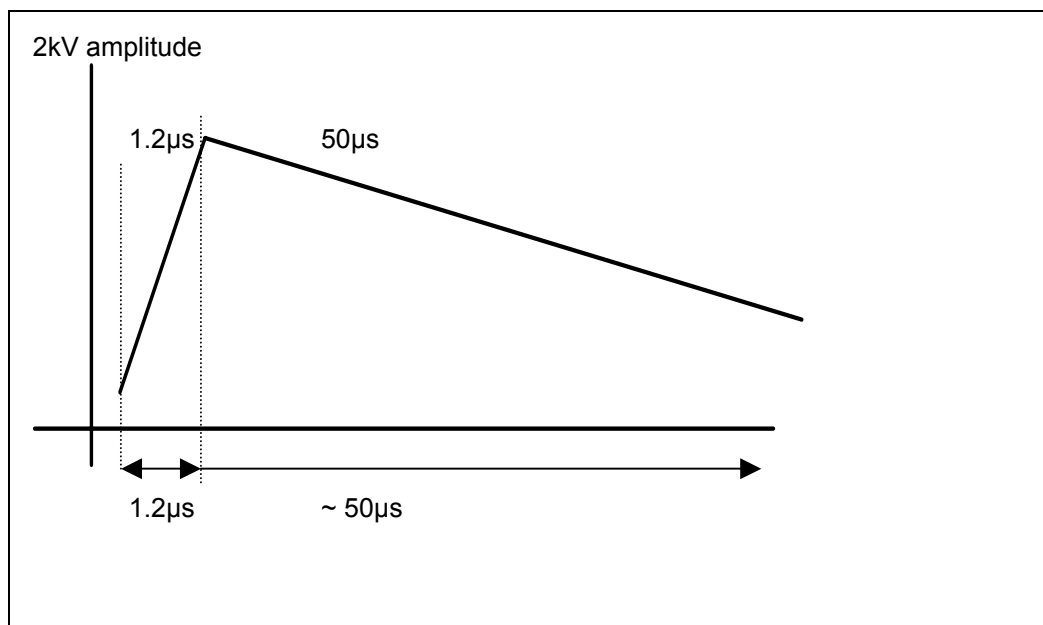

The Newsletter is for internal Philips use only. Issued by X-ray Generators XRD-Hamburg.
Hamburg, Germany. Tel: +49 40 5078 2084 Fax: +49 40 5078 1611
email: Guenther.Kramm@philips.com

XRD-Hamburg 074

CSIP Level 1 (05.0)

1/2

This document and the information contained in it is proprietary and confidential information of Philips Medical Systems ("Philips") and may not be reproduced, copied in whole or in part, adapted, modified, disclosed to others, or disseminated without the prior written permission of the Philips Legal Department.

(c) 2005 Koninklijke Philips Electronics N.V.
ALL RIGHTS RESERVED

The current of a varistor is 1 mA at rated voltage, with increasing voltage the current of the varistor rises according to the e-function.

Example of type S20K60 nominal data:

$V_{RMS} = 60V$ $V_{DC} = 85V$

Current 1mA ~ 60-70V, ~40mA at 120V.

In case there is an error which causes a current > 20mA in the varistor (caused by the high voltage between N and PE) the FI switch must respond immediately to interrupt the phases and N. In this way rotor control unit and generator are protected against destruction.

Failures of the rotor control unit are caused by mains conditions which are insufficient. If there are such problems with regard to N-PE it is necessary to install a mains adaptation transformer.

This transformer makes sure that the necessary mains conditions for the internal supply of the generator and the supply of external consumers are properly provided.

Please order the certain transformer for:

OPTIMUS RAD 9890 000 0230x

OPTIMUS RF/C 9890 000 0260x

CSIP Level 1 & Proprietary Notice:

This document and the information contained in it is strictly reserved for current Philips Medical Systems ("Philips") personnel, Philips licensed representatives and Philips customers who have purchased a valid service agreement for use by the customer's designated in-house service employee on equipment located at the customer's designated site. Use of this document by unauthorized persons is strictly prohibited. This document must be returned to Philips when the user is no longer licensed and in any event upon Philips' first written request.

Liability / Warranty Disclaimer:

Philips provides this document without warranty of any kind, implied or expressed, including, but not limited to, the implied warranties of merchantability and fitness for a particular purpose. Philips has taken care to ensure the accuracy of this document. However, Philips assumes no liability for errors or omissions and reserves the right to make changes without further notice to any products herein to improve reliability, function, or design. Philips may make improvements or changes in the product(s) or program(s) described in this document at any time.